

# HORNET

FIXED FRAME DAY WHEELCHAIR



The Melrose Hornet is available in Aluminium or Titanium and has been designed and built to be one of our lightest wheelchair's in the Melrose lifestyle range.

Primarily designed for active users, it is now also used by low

strength users for ease of pushing and lifting into vehicles.

The chair is completely custom built with camber, wheel sizes, seat width, depth and heights to your individual requirements.

Options include paediatric, and suspension versions.



# HORNET

## FIXED FRAME DAY WHEELCHAIR

Weight Limit: 100kg max, Custom Weight option also available  
Approx Weight: 4.6kg - Titanium standard 16" frame without wheels

### STANDARD SPECIFICATIONS

Titanium Frame 90460

Aluminium Frame 90460A

Adjustable Back Height and Angle

Fixed or Folding Back

Seat Width between 13" - 20"

Seat Length between 13" - 20"

Adjustable Foot-plate Height

Fixed Seat Height

Adjustable Balance Point

Brakes: Push to Lock Brakes  
Alloy Scissor Brakes

Castors: Roller-blade  
4", 5", 6" Solid  
4", 5", 6" Wide Solid

Quick release 24" Spoked Wheels on Alloy Rims

Tyres: 24 x 1 High Pressure Tread or Smooth  
24 x 1 Solid

Black Fabric Sling

Black Fabric Back with Leather Trim

Leg Strap

Five Year Frame Warranty

### OPTIONS AND ACCESSORIES

Castors: 4", 5", 6" Solid Alloy  
4", 5", 6" Soft Roll  
6", 8" Pneumatic

Castor Forks Frog Leg Suspension Forks  
Hybrid Frog Leg Suspension Forks  
Ultra Lite Forks with Bullet Pivots

Brakes: Composite Scissor Brakes  
Eagle Scissor Brakes  
Ds-Locks Brakes  
Brake Extension Levers  
Hill Climber Brakes

Wheels: Spoked 22", 24", 25", 26"  
BMX Wheels  
Xcore Wheels  
Spinergy Wheels  
Soft Wheels  
Loop Wheels

Push-rims: Rubber Coated - Red or Black  
Vinyl  
Foam  
CarboLife  
Natural Fit

Push handles: In Back Canes  
Removable  
Fold Down

Skirt Guards: Fabric Skirt Guards  
Fabric covered ABS Skirt Guards  
ABS Drop In Skirt Guards  
Alloy or Carbon Fibre Skirt Guards

Armrest: Tubular Swing Away  
Desk or Full Length Pad  
T Adjustable Height, with/without Skirt Guards

Anti-Tips with Plastic Wheels or Roller-blades

Adjustable Back Upholstery


Rear Wheel Camber 0° - 5°

Rear Suspension

Under Chair Bag

Fabric or ABS Foot-plate Cover

Cushion, Impact Guards, Castor Covers

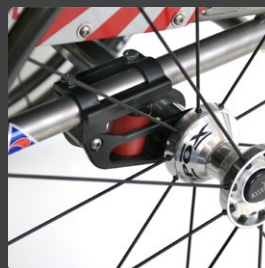

*The Hornet has suspension options for the rear wheels (left) - add 500gm  
Frog Legs Forks (centre) - add 220gm or Hybrid Frog Legs Forks (right) - add 160gm.*